



尺寸:65*140mm
 字体:思源黑体
 内容颜色:100%黑

4. Switch Mode

- ON/OFF**
Short press the key to start the toothbrush. After the toothbrush works in the selected mode for 8 seconds, and then short press the key to stop brushing.
- Function switching**
Within 8 seconds after the toothbrush starts in the selected mode, you can press the key to switch the different gears.
- Functional order**
Shortly pressing button switch can be switched in sequence: Gum Care,White,Clean.
- Booting cushioning**
1.5 seconds after power on, the intensity gradually increases to the standard intensity.

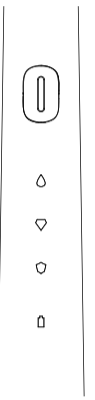

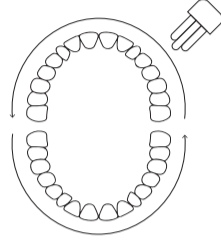


The function of automatic timing change zone reminder

- > 2 minutes timer: The dentist recommends brushing your teeth for 2 minutes at a time. Shut down after 2 minutes in any mode.
- > 30 seconds zone change reminder: Remind every 30 seconds during use in any mode. This is to remind users to change their modes and brush their teeth to ensure a balanced cleaning of the four modes of the mouth.

Automatic memory function
Toothbrush starts, after working in the selected mode for 8 seconds, short press the key to turn off the toothbrush, when you restart the toothbrush, it will automatically return to the mode before closing.

USER MANUAL

SN301

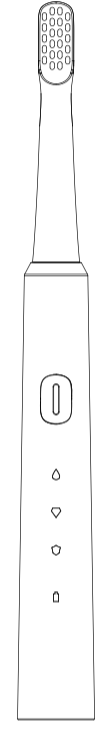


5. Instructions

Brushing steps


To prevent toothpaste from splashing, take a small amount of toothpaste on the brush head, soak it in water and then brush your teeth.
 Hold the toothbrush body gently, move slowly in a circular motion, do not press or scrub hard.
 Please rinse your mouth with water after brushing.

6. Charging instructions

The charging time is about 3.5 hours, and it can be used for two months after being fully charged (according to the cleaning mode, two minutes in the morning and night).

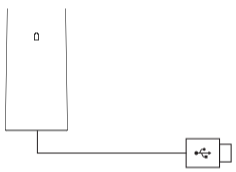

1. Product Description

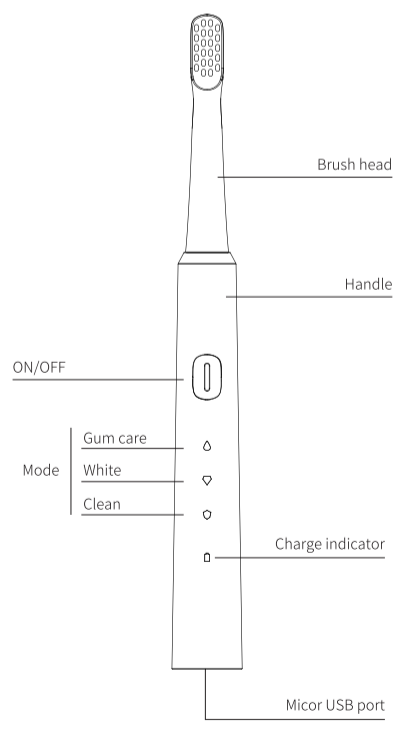

2. Specification

Model: SN301
 Demesion: 24.5X24.5X231.5mm
 Net weight: 60g
 Material: ABS+TPE
 Mode: Gum Care,White,Clean
 Waterproof: IPX7
 Input: Rated voltage 5V = Rated power: 5W
 Battery capacity: 750mAh
 Charger: Micro USB

Power Indicator


- > Regular using: white light is on
- > While charging: red light is on
- > When fully charged: grenn light is on
- > When battery is low: Red light flashes

Functional mode lights
In the power-on state, the corresponding mode light is solid white.
In the power-on state, the mode light corresponding to the indicator light is solid white.

Low battery protection released
After the low battery automatically shuts down, pressing the key again will not enter the working state. The battery indicator will turn into a red breathing light.

3. Install brush head

Hold the handle, insert the brush head 45 degrees to the right, and then rotate it to lock the left.

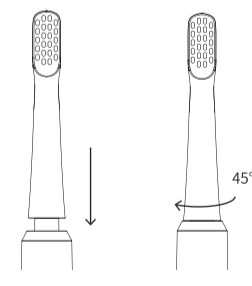

Replacement of bursh head
It's recommended to replace the toothbrush head every 2-3 months. When removing the toothbrush head, rotate 45 degrees to the right.

7. Maintenance

- Turn off the power supply of the electric toothbrush, rinse the handle and the brush head with water, tap the toothbrush head to ensure that there are no debris in the bristles, and dry the water on it.
- Wash the bristles and brush head after we use, do not pull the bristles to easily deform and cause shedding.
- Wash the toothbrush body with water, dry it with natural wind, do not rinse with water or acidic cleaning solution higher than 40 degrees.
- Dentist recommends replacing the brush head every 2-3 months.

1 Package: 1X handle, 1Xbrush head,1Xusb cable